




DNA Ladder 100 bp plus NO DYE

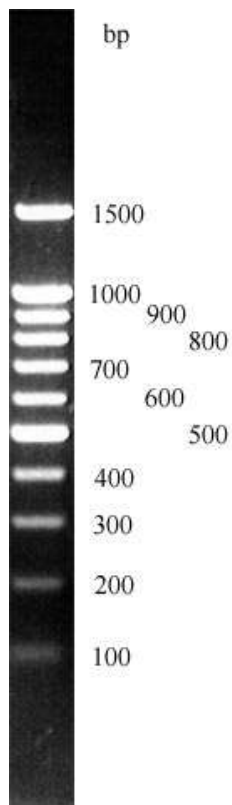

Features: GeneON 100 bp DNA PLUS Marker is composed of 11 individual DNA fragments which are re-dissolved in storage buffer.

Description/Preparation: The 100 bp PLUS DNA Ladder has a number of proprietary plasmids which are digested to completion with appropriate restriction enzymes to yield 10 bands suitable for use as molecular weight standards for agarose gel electrophoresis. The 100 bp PLUS DNA Ladder was designed for precise qualitative and approximately quantification of DNA mass. There are no unspecific bands besides the fragments. The 500bp fragment is used a reference band.

Usage: 1 µg/lane

Concentration: 0,2 µg/µl

Storage Buffer: 10 mM Tris-HCl (pH 8.0), 1.0 mM EDTA, 50mM NaCl

Tracking Dye: 6X DNA Loading dye required

Number of bands: 11 fragments

Gel-picture: 1 µg of 100bp + 1,5Kb DNA Ladder visualized by ethidium bromide staining on 1,7 % TAE agarose gel

Number of bands and approx. DNA Mass:

base pairs	100	200	300	400	500	600	700	800	900	1000	1500
DNA mass ng	20	20	30	50	150	80	70	110	120	200	150

Loading: 0.5 - max. 1.2 µg/ lane in 1.7 - 2.5 % Agarose. We recommend to add the loading dye in amount 1/10 from the DNA Ladder volume

Transportation: shipped on blue ice or room temperature

Storage: at -20°C for 24 months

Ordering information:

Cat.-no	Description	Amount
306-005	100 bp DNA Ladder PLUS, no stain	50 µg/250 µl
306-025	100 bp DNA Ladder PLUS, no stein	5 x 50 µg

DNA Ladder 100 bp plus NO Dye

GeneON .. a good decision ..

Contact Phone +49-(0)-621- 5720 864 Fax: +49-(0)-621-5724 462


E-Mail: <mailto:info@geneon.net> WEB: <http://www.taq-dna.com/> Version: 25.10.2009 AS

Unless specified otherwise, all products of GeneON are sold for research use only.

.. a good decision ..



Cat. No. 306-005, 50 µg

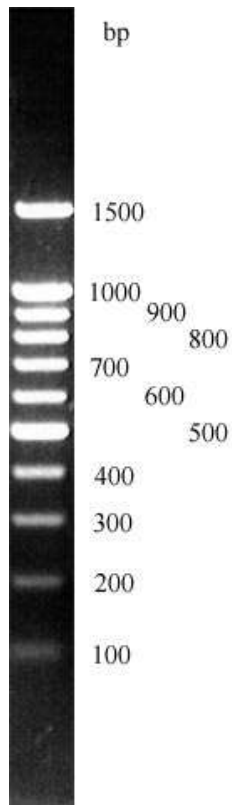

Features: GeneON 100 bp DNA PLUS Marker is composed of 11 individual DNA fragments which are re-dissolved in storage buffer.

Description/Preparation: The 100 bp PLUS DNA Ladder has a number of proprietary plasmids which are digested to completion with appropriate restriction enzymes to yield 10 bands suitable for use as molecular weight standards for agarose gel electrophoresis. The 100 bp PLUS DNA Ladder was designed for precise qualitative and approximately quantification of DNA mass. There are no unspecific bands besides the fragments. The 500bp fragment is used a reference band.

Usage: 1 µg/lane

Concentration: 0,2 µg/µl

Storage Buffer: 10 mM Tris-HCl (pH 8.0), 1.0 mM EDTA, 50mM NaCl

Tracking Dye: 6X DNA Loading dye required

Number of bands: 11 fragments

Gel-picture: 1 µg of 100bp + 1,5Kb DNA Ladder visualized by ethidium bromide staining on 1,7 % TAE agarose gel

Number of bands and approx. DNA Mass:

base pairs	100	200	300	400	500	600	700	800	900	1000	1500
DNA mass ng	20	20	30	50	150	80	70	110	120	200	150

Loading: 0.5 - max. 1.2 µg/ lane in 1.7 - 2.5 % Agarose. We recommend to add the loading dye in amount 1/10 from the DNA Ladder volume

Transportation: shipped on blue ice or room temperature

Storage: at -20°C for 24 months

Ordering information:

Cat.-no	Description	Amount
306-005	100 bp DNA Ladder PLUS, no stain	50 µg/250 µl
306-025	100 bp DNA Ladder PLUS, no stein	5 x 50 µg

.. a good decision ..